




## Instructions for the Baxter Mini-Infuser™ Pump 300XL

### Before you begin:

1. Wash your hands.
2. Gather your supplies:
  - 2 Normal Saline syringes (remove air)      □ Heparin syringe
  - Alcohol pads
  - Baxter Mini-Infuser Pump
  - Your medication syringe


### Giving Your Medication:

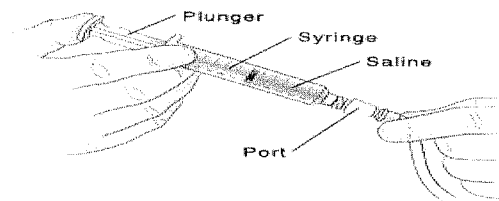

#### Flush your access device

- Disinfect the injection cap attached to your access device by holding the end of the injection cap with your alcohol wipe and using a twisting motion to scrub. Scrub for 15 seconds. Allow to dry.
- Open any clamps on access device tubing, remove the protective cap on the syringe and flush with \_\_\_\_\_ cc of Normal Saline
- Discard the Normal Saline syringe


#### Prime your tubing

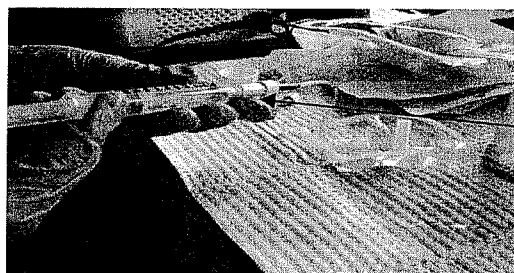


Syringe Hub

- Attach infusion tubing to **syringe hub**
- Push **syringe plunger** gently until the tubing fills completely with medication
- Set aside in safe location

## Administering Your Medication:

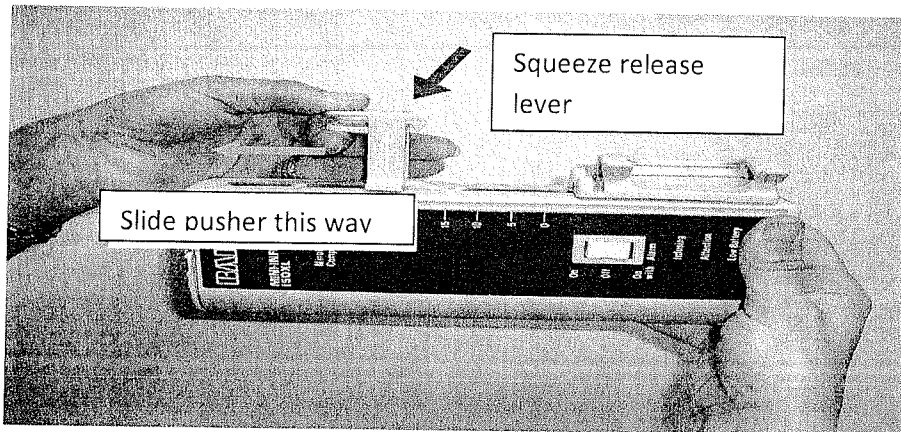

1. Grasp **Release Lever** on the Baxter Mini-Infuser pump. Squeeze in direction of arrow until the **Pusher** moves freely, slide to top of pump.

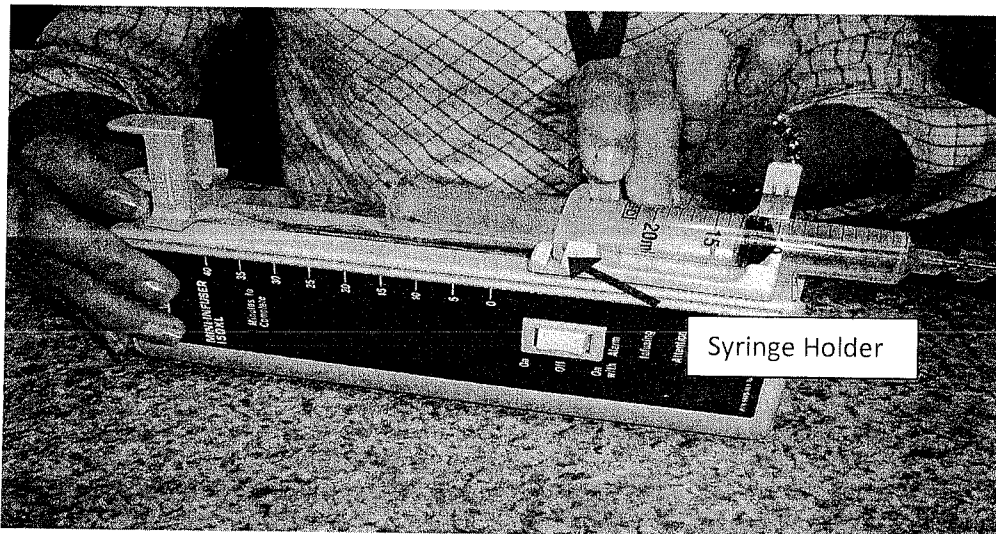

2. Snap syringe barrel into the **Syringe Holder** of the Baxter Mini-Infuser. Grasp **release lever** and slide **pusher** until it grasps the top of the **syringe plunger**, then release the lever.
3. Place the pump in a convenient location.
4. Attach the primed tubing to your access device
5. Start your infusion by turning the **On-Off/Switch** to the "On with Alarm" position. You will see a green light flashing when infusing correctly.

6. When your infusion is complete turn the **On-Off/Switch** to the "Off" position. Squeeze the **release lever** completely and reset the **plunger** to the top of the pump. Remove the empty syringe and dispose.